



NATIONAL SURGICAL CORPORATION

CSSD Height Adjustable Packing Tables

ISO 9001
BUREAU VERITAS
Certification



A proudly founded and owned Australian company - providing value to our customers beyond the product

CSSD Height Adjustable Packing Tables



National Surgical Corporation Height Adjustable Tables: high quality, customised, adjust height as required

AUSTRALIAN DESIGNED & MANUFACTURED

Designed to meet your individual needs and manufactured locally in Australia

VARIOUS SIZES AVAILABLE

Available in various sizes up to 2.3 metres wide. Comes standard with castors and electric linak lifting legs

HIGH QUALITY STAINLESS STEEL

Manufactured with 304 grade stainless steel ensures a high quality and durable finish

IMPROVE PRODUCTIVITY WITHIN CSSD

Improve productivity with your sterilisation department by allowing the table to electrically adjust the benchtop height to suit all users. Range of accessories improve workplace functionality

REDUCE POTENTIAL OH&S RISKS

Press button control to raise or lower benchtop to required height (750mm - 1400mm) required by the user, to reduce the risk of staff OH&S issues

IDEAL FOR CSSD

Designed especially for any hospital, clinic or day surgery sterilisation/CSSD department

SAFE WORKING LOAD


Each table is able to withstand up to 100kg

EXPERTS IN CUSTOM DESIGN

NSC has designed and manufactured a range of shelving & storage solutions for CSSDs including wrap racks, wire racking, chrome wire shelving and compactors



LINAK LIFTING LEGS

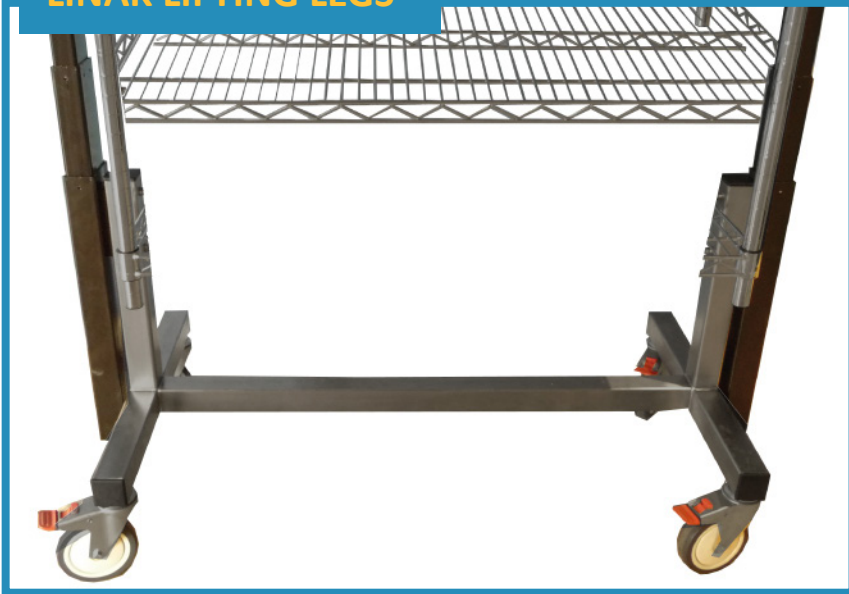

The high quality LINAK® electric lifting legs allows for a seamless and low noise lifting of the stainless steel table top.

The LINAK® lifting legs are mounted to sturdy side columns which ensures the durability and strength of the table.


The LINAK® legs feature an easy to use up/down arrow control panel (mounted to the front side of the table).



HEAVY DUTY CASTORS

The NSC CSSD height adjustable packing tables are manufactured with heavy duty castors to allow for support and weight loading of the table.

These high quality castors come standard as 2 swivel and 2 braked to lock the table in place. Tables can also be manufactured as static height adjustable tables with feet.

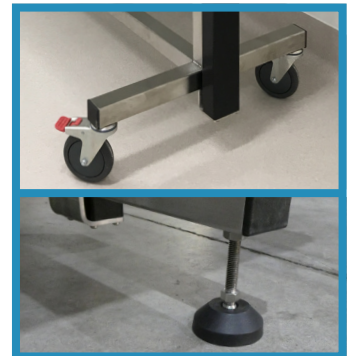

DIMENSIONS

The NSC CSSD height adjustable packing tables are custom made to order.

Tables can be manufactured in any size - from as small as 550x750mm to as large as 2400x900mm.

Below are the most common sizes:

1200 x 750mm	2000 x 750mm
1200 x 900mm	2000 x 900mm
1500 x 750mm	2300 x 750mm
1500 x 900mm	2300 x 900mm
1800 x 750mm	
1800 x 900mm	

Packing tables can come in any size to suit your sterilisation department or packing area



CSSD Height Adjustable Packing Tables

ACCESSORIES

UNDERBENCH SHELF

Recommended for all tables, the underbench chrome wire shelf and rear overbench shelf allows for storage of peel packs, consumables and equipment when not in use. Size of shelf is dependent on size of table - can be single or double shelves.



Table size	Recommended underbench shelf size
1200 x 750mm	910 x 530mm
1200 x 900mm	910 x 610mm
1500 x 750mm	1070 x 530mm
1500 x 900mm	1070 x 610mm
1800 x 750mm	1360 x 530mm
1800 x 900mm	1360 x 610mm
2000 x 750mm	1520 x 530mm
2000 x 900mm	1520 x 610mm
2300 x 750mm	1820 x 530mm
2300 x 900mm	1820 x 610mm

OVERBENCH REAR SHELF

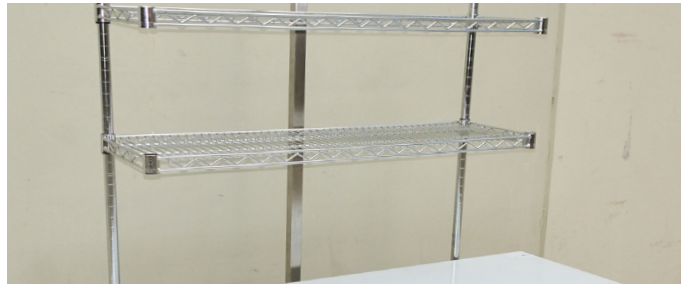


Table size	Recommended overbench shelf size
1200 x 750mm	1070 x 300mm
1200 x 900mm	1070 x 300mm
1500 x 750mm	1220 x 300mm
1500 x 900mm	1220 x 300mm
1800 x 750mm	1520 x 300mm
1800 x 900mm	1520 x 300mm
2000 x 750mm	1520 x 300mm
2000 x 900mm	1520 x 300mm
2300 x 750mm	1520 x 300mm
2300 x 900mm	1520 x 300mm

MONITOR ARM HOLDER

Allows for screens and monitors up to 32" to be mounted to the table at eye level. Easy one-hand repositioning with an inbuilt mechanical spring. It frees up table space and enhances ergonomics.

Weight loading: 8kg max



Code: HAT-C DS112D
Monitor arm holder

KEYBOARD HOLDER

Combine the monitor arm holder with a retractable keyboard holder for a complete technology user experience.

Store keyboard under table and retract to allow for quick and easy typing. Retractable mount allows any sized keyboard to be stored under table when not in use. Weight loading: 2kg



Code: C KEYBOARD DR
Keyboard holder retractable with mount - 660x343x1.5mm



PUREKEYS FLAT MEDICAL KEYBOARDS

- Completely flat silicone top
- Easy to clean and disinfect with all hospital grade disinfectants
- Water and spillproof, no liquids can penetrate the keyboard
- Connects easily with USB wired cable
- Fixed angle option to allow keyboard to be permanently viewed and used on an angle
- Reduce contamination in the hospital, dust and debris will not get trapped, like in raised keyboard keys
- Place on table or use with retractable keyboard mount



Code	Description	Image
K-30000100	Medical keyboard compact 103 keys IP66, wired	
K-40000100	Medical keyboard compact 104 keys, fixed angle, IP66, wired	
K-40100100	Medical keyboard compact 87 keys with touchpad, fixed angle, IP66, wired	
K-50500100	Medical Keyboard full size 105 keys IP66, wired	
K-40065	Medical mouse, touch scroll, IP66, USB	
K-90005002	4 self sticking magnets and 4 metal discs for fixing 300 series	

Wireless versions are also available

CSSD Height Adjustable Packing Tables

MAGNIFIER LIGHTS

Magnifying LED lights allow sterilisation technicians to perform their work to the greatest accuracy when viewing instruments and equipment when being cleaned and packed. The LED generates light without any dark spots and the magnifier allows for enhanced inspection of the items.



Code: MEK DAYA25130

Daylight Omega 7 Magnifier Glass LED

- 175mm lightweight acrylic lens and high-power LEDs
- 4 brightness levels for clear viewing
- 3 dioptres (1.75X)
- Internal spring arm allows accurate positioning
- Long reach 1180mm
- LED rated life (hours): 50,000



Code: MIM-09340


Lupa LED H.F

- 127mm LED magnifier with round lens
- Adjustable height articulated arm
- 5 dioptres
- Magnification lens diameter 127mm
- Articulated arm length: 900mm
- LED rated life (hours): 40,000
- Various mounting options

ROLL HOLDER

Securely store steri-peel rolls, tapes and other rolls at rear of table to allow for easy access and cutting.

Available as single or double holders. Size depends on the size of the table and if there are rear shelves included

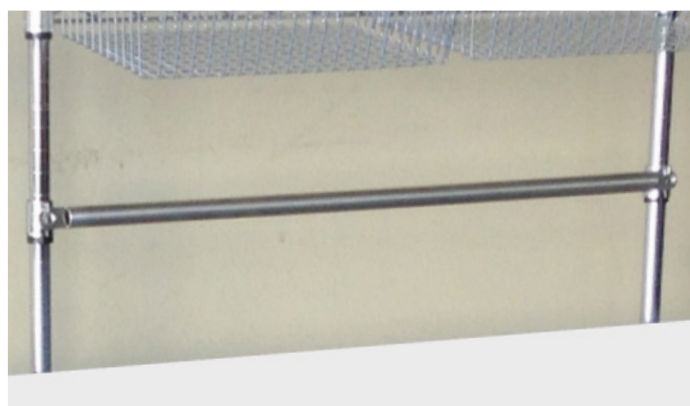

Code: C PAPER ROLL HOLDER


Removable roll holders with mounts on each end

Available sizes:

- 910mm
- 1070mm
- 1220mm
- 1360mm
- 1520mm
- 1820mm

CONVEYOR ROLLER

Improve OH&S of sterilisation personnel by moving heavy trays from table using the steel conveyor roller. Trays can be easily rolled onto trolleys or lifters. Conveyor rollers can be mounted on any side of the table

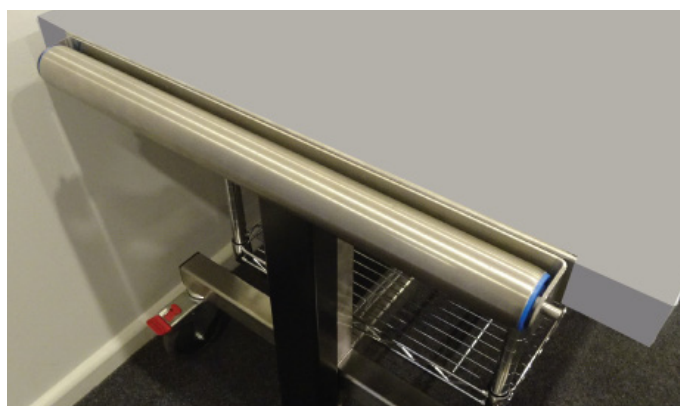

Code: C ROLLER

Available sizes:

- 290mm
- 450mm
- 600mm



UNDERSHELF LIGHT STRIP

LED light strip mounted under rear shelf provides superior brightness to the table. Turn off easily with close-by switch when not in use. Size of strip depends on rear shelves installed



Code: C LED LIGHT

POWERBOARD

When additional power requirements are needed for computers, monitors, keyboards, heat sealers and other devices. The powerboard is mounted under rear of table keeping cables away but still accessible.



Code: C POWERBOARD

Available with 4, 5, 6 power outlets

FAMOS ROLL HOLDERS & CUTTERS

Various roll holders with cutting device to create the length of packaging material required for instruments being sterilised. The rotary knife of the cutting device is made from the highest quality hardened steel and is self-sharpening which ensures years of maintenance-free use.

- Nylon bearings makes it easier to unroll the steripeel rolls
- Top extension option doubles storage capacity


Code	Description	Extension code	Image
FM-50.104.259 (TS450)	Roll holder with cutting device Storage capacity: 450mm dimensions: 450 x 408 x 112mm	N/A	
FM-50.100.140 (TS600)	Roll holder with cutting device Storage capacity: 600mm dimensions: 820 x 420 x 130mm	FM-50.100.138	 <i>pictured with extension</i>
FM-50.100.080 (TS750)	Roll holder with cutting device Storage capacity: 750mm dimensions: 970 x 420 x 130mm	FM-50.100.081	 <i>pictured with extension</i>
FM-50.100.130 (TS900)	Roll holder with cutting device Storage capacity: 900mm dimensions 1120 x 420 x 130mm	FM-50.100.131	
FM-50.101.118 (TS1200)	Roll holder with cutting device Storage capacity: 1200mm dimensions: 1200 x 330 x 170mm	FM-50.101.269	

CSSD Height Adjustable Packing Tables

FAMOS HEAT SEALERS

Redesigned and innovative PROTEC heat sealers by Famos. The housing of the Famos rotary heat sealers has been designed with a hygienic point of view. Surface disinfection and cleaning can be easily done due to the round corners and closed housing. The machines can be easily opened and serviced without the use any tools. Feed-in guide integrated into the machine can be quickly adjusted by an easy clicking system (0-40 mm).

F108 PROTEC series rotary heat sealers

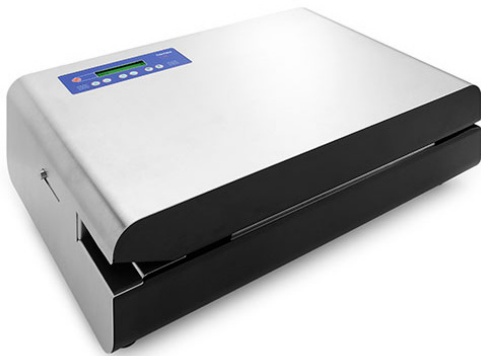

- Basic model
- No internal printer
- 3 models in range (distinguish by a menu controlled operation and micro processor controlled system)
- Digital temperature setting and control
- LCD display
- USB connection
- Automatic transport start and reverse transportation feature prevents instruments from getting damaged
- Safety feature stops the machine from sealing when there is a temperature difference of 5°C against set temperature

Code	Description
FM-30.100.550	Famos F108 PROTEC rotary heat sealer
FM-30.100.551	Famos F108 Tyvek PROTEC rotary heat sealer
FM-30.100.554	Famos F108 TX PROTEC rotary heat sealer
FM-30.100.555	Famos F108 TX Tyvek PROTEC rotary heat sealer
FM-30.100.558	Famos F108 TX-P/USB PROTEC rotary heat sealer
FM-30.100.559	Famos F108 TX-P/USB Tyvek PROTEC rotary heat sealer
FM-30.100.623	Famos front table for F108 PROTEC series
FM-30.100.620	Famos roller front table for F108 PROTEC series

F110 PROTEC series rotary heat sealers




- Parameter control has been integrated - continuously monitors the parameters temperature, pressure and speed to be able to validate sealing process
- If a parameter is incorrect, a warning message with sound alarm will be displayed and the machine will stop sealing
- Monitoring of the parameters can be done in real time on a PC
- 3 models
- Automatic transport start, reverse transportation feature and a production counter with alarm
- Integrated internal printer which can print various data on the pouches (packing date, expiry date and a barcode)
- Can be connected to various instrument tracking systems
- Available with a grooved pressure roller (rib seal) for laminated or paper pouches or a pressure roller without grooves (flat seal) which is especially suitable to seal Tyvek

Code	Description
FM-30.100.570	Famos F110 D PROTEC rotary heat sealer
FM-30.100.571	Famos F110 D Tyvek PROTEC rotary heat sealer
FM-30.100.574	Famos F110 DEP-P/USB PROTEC rotary heat sealer
FM-30.100.575	Famos F110 DEP-P/USB Tyvek PROTEC rotary heat sealer
FM-30.100.582	Famos F110 DBC-P/USB PROTEC rotary heat sealer
FM-30.100.583	Famos F110 DBC-P/USB Tyvek PROTEC rotary heat sealer
FM-30.100.626	Famos front table for F110 / F220 PROTEC series
FM-30.100.619	Famos roller front table for F110 / F220 PROTEC series

CSSD Height Adjustable Packing Tables


F220 PROTEC series rotary heat sealers

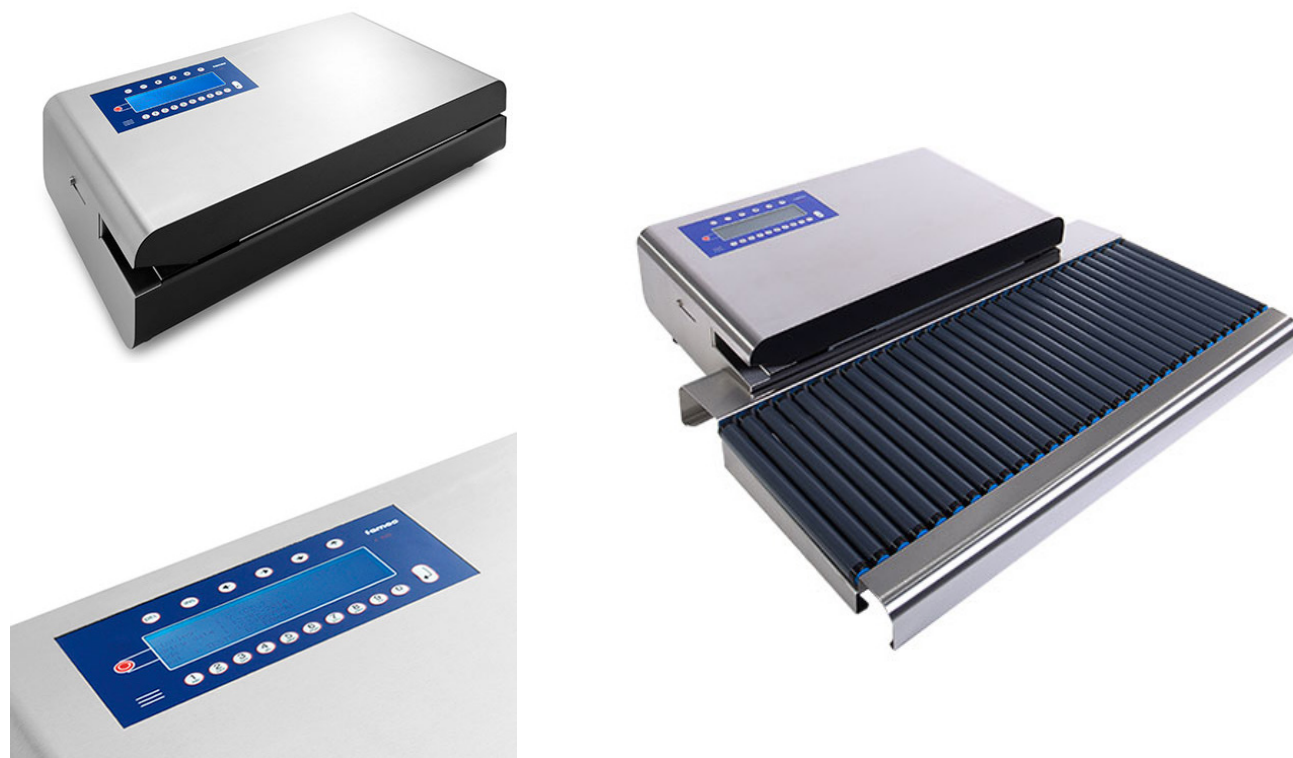

- Most extensive sealers of the Famos sealer range
- 2 models - equipped with parameter control which makes the sealers DIN ISO 11607-2 compliant
- Parameters temperature, speed and pressure are constantly being checked with each seal
- If one of these parameters is incorrect, the machine will stop sealing and will notify the user by a warning message and a sound alarm
- Software program is supplied
- Integrated internal printer which print the packing date, expiry date, a code and one line of text which can consist of up to 30 characters + barcode and tracking data

Code	Description
FM-30.100.600	Famos F220 DEP-P/USB PROTEC rotary heat sealer
FM-30.100.601	Famos F220 DEP-P/USB Tyvek PROTEC rotary heat sealer
FM-30.100.608	Famos F220 DBC-P/USB PROTEC rotary heat sealer
FM-30.100.609	Famos F220 DBC-P/USB Tyvek PROTEC rotary heat sealer
FM-30.100.626	Famos Front table for F110 / F220 PROTEC series
FM-30.100.619	Famos Roller front table for F110 / F220 PROTEC series

F70 series impulse sealers

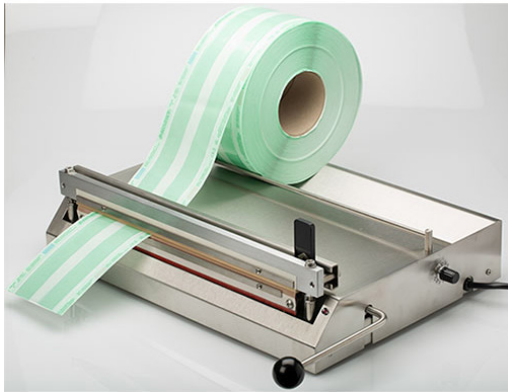

- Suitable for small quantities of seals
- Seal pouches of polypropylene, polyethylene, Tyvek, paper pouches with glue strip, dust covers and other materials
- Completely made with AISI 304 stainless steel
- Ideal for processing packaging material from the roll with its built-in roll holder and a cutting device
- Sealing time and cooling down time are being shown on a LED
- Time-controlled and produce a seal width of 8mm
- F70-300 has an effective seal length of 320mm
- F70-400 has an effective seal length of 420mm

Code	Description
FM-31.100.008	Famos F70 impulse sealer, 300mm time control
FM-31.100.016	Famos F70 impulse sealer, 400mm time control


RACK/SUPPORT TABLE

Stainless steel table - 250x175x45mm.

Stainless steel rack to raise equipment off the table top. Use as a support ledge or to raise monitors or other electrical equipment to protect it from water which may be produced during the cleaning process.

Other sizes available.

Code: TBL-SS-100

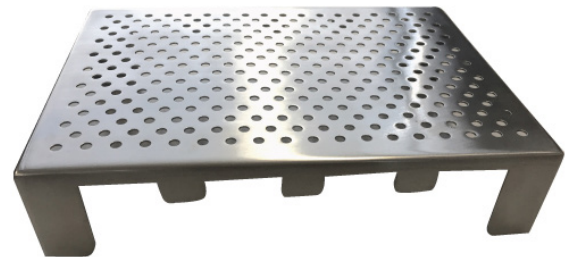

INSULATION TESTERS

- Compact, handheld, battery operated insulation tester
- Used to detect and locate defects such as pinholes, cracks and bare spots in protective coatings applied over conductive instrument's surfaces
- Insulation defect testing of electrosurgical instruments and cables



Code	Description
LTA-MM513-600	Electrosurgical insulation defect detector
LTA-MMSBT-170	Saddle block
LTA-MMTRI-0022A	Tri hole electrode
LTA-MMLSE-0029	LS ring Electrode
LTA-MMBRU-0007	Brush Electrode
LTA-MMGWC-0005L	Green ground wire with clamp
LTA-MM513-120EU	AC charger
LTA-MMWIT-200A	Wire tester
LTA-MMBPT-190	Bi-polar fixture



Let NSC Support You

**BE SUPPORTED BY AN
AUSTRALIAN COMPANY**

PROVIDING VALUE BEYOND THE PRODUCT

What makes NSC different to our competitors to support you?

- A proudly Australian owned and founded company
- Combining 25 years of expertise in the healthcare industry
- Creating value for our customers beyond the product
- Working with our customers to offer specialised solutions
- Offering a large product range focused on cost effectiveness, patient care and productivity
- Low overheads allows us to provide cost savings to our customers compared to the multinationals
- Coordinated, whole company approach to every customer
- Establishing trust for long term customer relationships

Who are we?

National Surgical Corporation (NSC) is an import manufacturer in the healthcare industry, distributing surgical instruments, consumables, equipment and shelving to all major hospitals and healthcare facilities throughout Australia. Our extensive product range can be used by all specialist disciplines and departments within a hospital, day surgery or specialist centre.

Let NSC support you and your department

Don't be afraid to change, NSC will support and assist you throughout every step. NSC has an extensive training and inservicing material available and can tailor a program to meet your individual needs so that all personnel are aware of how the product operates. Switching to NSC can save you significantly and we can show you how!

Our range includes:



Wound Drains




Electrosurgical Range



OT & General Consumables



Surgical Instruments

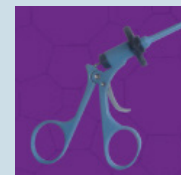

Laparoscopic Instruments




Medichill Ice & Heat Packs



Protection Products



Cleaning & Sterilisation

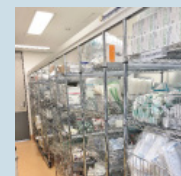

Shelving & Storage



Pegasus High Density Storage

CUSTOMER SERVICE

Contact your National Surgical Corporation Account Manager for more information

Phone: 1800 138 138 (free call)

info@nationalsurgical.com.au

www.nationalsurgical.com.au